

DENTAL BIB



Our **Dental Bib** is designed with unique bonding technique which bonds the layers of tissues and poly securely together, eliminating separation of the layers. The horizontal embossing allows the moisture to flow across the towel and securely locks the liquid deep at bottom.

Made of 2 ply paper and 1 ply PE. Extra heavy weight 18+18gsm origin paper provides excellent absorbency.

Size of 45cm x 33cm

Code : 5-800
Color : Blue, green, lavender, yellow, pink
Packing : 50pcs x 10packs



CROSS PROTECTION®



DENTAL & MEDICAL DISPOSABLE



CROSS PROTECTION (M) SDN BHD
No.11, Jalan Perindustrian 5, Batu 5 1/2,
Off Jalan Haji Abdul Manan, 41050 Klang,
Selangor, Malaysia.


✉ operation@crossprotection.com
marketing@crossprotection.com
sales@crossprotection.com

☎ +603 3392 0000 / +6016 262 0436

Made In Malaysia

www.crossprotection.com.my

SALIVA EJECTOR

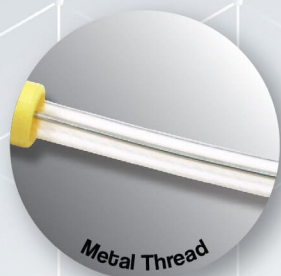

Metal Thread



Smooth Edge



Bend Easily • Stay & right at the position



Smooth Suction • No Block

Our Saliva Ejector 145mm is made of transparent, atoxic PVC. The selected raw materials, patented design featuring round edges, smooth surfaces and perfectly placed suction ports that help to enhance patient comfort while protecting delicate mucosal tissue.

INTENDED USE:

The Cross Protection® SALIVA EJECTOR are intended to be used as suction tubing of removing fluids and saliva from the open mouth during dental procedure without causing any abrasion or tear to the soft tissue.

INSTRUCTIONS FOR USE:

- 1) Take the saliva ejector, verify the safety and the functionality of the device, and assure yourself that the ejectors are integral; if the ejector has a fixed tip verify that it is effectively fixed.
- 2) Bend the saliva ejector to your specific needs and pay attention to not crush too much of the tube and check if something undesired happened (discharge of the inner wire, loss of the tip or something else).
- 3) Put the saliva ejectors in the specific place on the tube of the aspirator.
- 4) Put the saliva ejectors into the mouth of the patient and make sure to use it with care and diligence.
- 5) If you are using a saliva ejector with a removable tip pay attention to the accidental detachment of it during the various operations inside the mouth of the patient. The images given are only examples and they are not exhaustive.

SIZE	PACKING	CODE	TIP COLOR	TUBE COLOR
145mm	100pcs x 10bags	5-501MX	Assorted Color	Clear
		5-501BB	Blue	Blue
		5-501BC	Blue	Clear
		5-501W	White	White
		5-501CW	White	Clear
		5-501GG	Green	Green
		5-501PP	Pink	Pink
		5-501Y	Yellow	Yellow

MEDICAL DISPOSABLE

Why use disposable?

Disposable clothing, linens and other items dramatically decreases time consuming chores. No more time or money wasted on expensive washing, bleaching, drying or ironing procedures. You can also say goodbye to multiple washings of the same load because some items remain unwashed. Our linens, drapes and clothing always come new, hygienic, clean, relatively cheap and available for immediate use.

Does disposable save money?

We sell products that cost roughly the same amount or even less than the cost of washing similar, traditional materials. When one includes the other products and materials that must be bought for the washing procedure, such as machine costs, fabric softener, water and electricity. You do not need to worry about saving money by purchasing of disposable cloth. To control your budget, you must know your various expenses and when they have to be paid. Disposable products can help you control your budget because they become a regular, measurable cost. You no longer have to fear staff falling sick or absent, unable to wash the linens. Disposable is always available and ready, at a lower cost than most reusable products when one factors in their care and cleaning.

Does disposable look professional?

Don't forget that image means a lot when it comes to clients. Always clean and hypoallergenic, disposable linens and clothing give an impression of impeccable hygiene. With disposable bedding, all your sheets are new and spotless from their first moment of use. Disposable makes your team more flexible and ready to respond to rapid increases in activity without worry. It's a form of thinking ahead.

Why choose Cross Protection?

We manufactured all kinds of linens, and apparels. Our trademark material Cotton Soft guarantee unsurpassed, overall performances in the areas of skin dryness, odor reduction and inhibition of bacterial growth.




CONSUMER
Uninterrupted sleep & reduced health risks for wounds and falls

CAREGIVERS
Improved morale & productivity



FACILITIES
Reduced product usage, laundry & odors

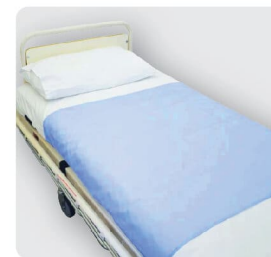

AAMI BLANKET

Code : 2-502
Size : 141x152cm
Color : Blue
Packing : 100pcs / carton



AAMI BEDSHEET

Code : 2-501
Size : 141x220cm
Color : Blue
Packing : 1pc x 100 bags



AAMI PILLOW CASE

Code : 2-504
Size : 74x48cm
Color : Blue
Packing : 10pcs X 20bags



COATED BEDSHEET

Code : 2-503
Size : 141x220cm
Color : White
Packing : 1pc X 100bags



COATED PILLOW CASE

Code : 2-500
Size : 74x48cm
Color : White
Packing : 10pcs X 20bags

Features:

Bacteria - Free*

The promise of unsurpassed overall performance in the areas of skin dryness, odor reduction, urine pH neutralization and inhibition of bacterial growth.

Fitted Corner *

Secures product in place (Applicable for Bedsheet only)

Cloth-Like Back Sheet *

Product outer layer is soft and rustle free.

Breathable*

Allows air to circulate around the body, promoting skin integrity.

Latex-Free*